

Basic Ankle Joint



Features

Small size for **toddlers**

+ Maintenance free, sustained joint friction **resists foot drop** during swing phase without springs

+ **Adjustable friction** for a smooth heel rocker

Fabricated from **high-strength aluminum**

Quick changes to alignment and movement stops with self-locking setscrews

Range of motion up to **20° dorsiflexion** and **20° plantarflexion**

Quiet movement stops promote acceptance

Lightweight aluminum stirrup

Standard

Fully assembled as one item for clear cost estimates and uncomplicated ordering

Integrated upright bar ensures **optimal transfer of force**

Upper section **CNC-machined from one piece** of material for increased stability

Precise assembly calibrated to 1/100mm with self-lubricating washers for smooth operation free of play

Maintenance-free stainless steel joint assembly with bronze bushing

Unrivalled **20-month guarantee**

Engraved **serial number** for traceability

Clearly illustrated instructions for use


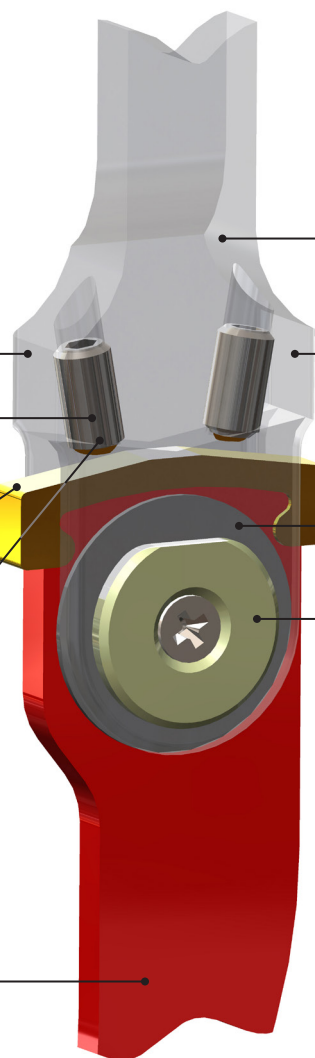









Fig. left

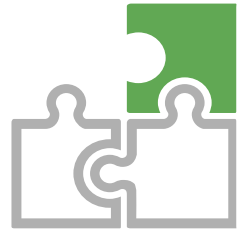

We recommend using our configurator to make your selection.
The configurator considers, among other factors, weight to height ratios.



Weight Classification (body weight)			Joint	
Article No.	AFO/KAFO	HKAFO/RGO	Width	Thickness
Basic-18	up to 10 kg	up to 10 kg	18 mm	8,0 mm

 <p>Adjustable Alignment</p>	Modular	Force Options	 <p>Range of Motion</p>	 <p>Swing Phase</p>	 <p>Quick Adjustment</p>	Interchangeable Stops	 <p>No Grinding Necessary</p>	 <p>Noise Dampening</p>	 <p>Water Resistant</p>
---	---------	---------------	--	---	---	-----------------------	--	--	--